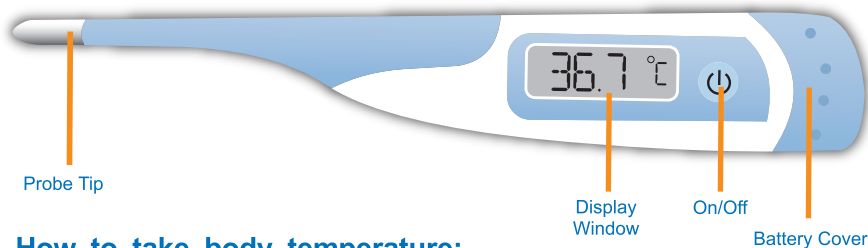






ThermoFlex™-10



How to take body temperature:

- 1 Pull out the thermometer from its protective cover.
- 2 Press the  button then the thermometer will “beep” and display all signals.
- 3 Next, it will display the last temperature reading for about 2 seconds. After that, temperature is being measured while  is flashing. Measuring time takes about 10 seconds.
- 4 When the highest temperature is reached and stabilized,  stops flashing and Bi — Bi — Bi — Bi — sound is heard. (short Bi-Bi-Bi — Bi-Bi-Bi — Bi-Bi-Bi is heard when temperature exceeds 37.8°C(100°F))
- 5 Remove the thermometer and take the reading. If the thermometer is not removed from the measured place (such as mouth, rectal or underarm) after the Bi sound, the indicated temperature will fluctuate slightly if a higher temperature is detected.
- 6 Press the  button to turn off the power or it will turn off automatically after 1 minute. Put the thermometer back into its protective cover.

For Oral Use – Place the probe tip well under the tongue and keep the mouth closed.

For Rectal Use – Lubricate the probe tip with lubricating jelly for easier insertion and gently insert the tip about 1/2” into the rectum. A rectal temperature is 0.5°C (0.9°F) higher than an oral temperature.

For Armpit Use – Wipe underarm with a dry towel, place the probe tip in armpit and hold firmly. An armpit temperature is 0.5°C (0.9°F) lower than an oral temperature.

Safety information:

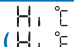


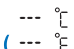
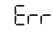
- High, prolonged fever requires medical attention.
- Do not allow children to take their temperature unattended. Stay still while taking temperature. The probe could present a hazard.
- Keep the thermometer away from children. Loose parts swallowed would be harmful.

ThermoFlex™-10




Useful tips:

- What is “Normal Temperature”?
37°C or 98.6 °F is generally accepted as “normal temperature”. However, 36.1°C (97°F) to 37.8°C (100°F) may still be considered as normal according to each individual’s age and average temperature while they are well.
- Do not drink hot or cold fluids, smoke or exercise prior to taking reading.
- Do not drop, bite or bend the probe beyond 45°.
- Keep the thermometer in its protective cover while not in use. Store in cool dry place.
- For personal hygiene, wipe / clean the thermometer with 70% Isopropyl alcohol before and after use.

Error Indications:

Error	Problems	Solutions
	Temperature $\geq 43^{\circ}\text{C}$ (109.4°F),	out of measuring range, incapable to take measurement.
	Temperature $\leq 32^{\circ}\text{C}$ (89.6°F),	out of measuring range, incapable to take measurement.
	Low battery indication,	replace the battery immediately.
	Incomplete measurement - If you take the thermometer out from measured place or move before beep sound is heard, “---°C(°F)” will come out with beep sound,	turn off the thermometer and then turn it on again to take a new measurement.
	Problem exists and can not take measurement,	send back to the local distributor for check up and repair.

Switch between °C and °F

Turn off the power and keep pressing the  button for 5 seconds. Then the LCD window will display the last reading in current scale. When C(F) appears on the screen, release the button and choose C or F by pressing the  button again. When the scale has been chosen, wait until  flashes to take a new measurement.

Replacing the battery—when “” appears on the LCD window:

- 1 Pull out the battery cover.
- 2 Remove the battery carefully from its holder by using a clip or a pointed object.
- 3 Insert a new battery (LR41) with the “+” symbol face down.
- 4 Put the battery cover back in place.

** Do not recharge the battery.

** Do not swallow the battery which may be fatal.

** Do not attempt to disassemble the unit except battery replacement.